
# ecommed

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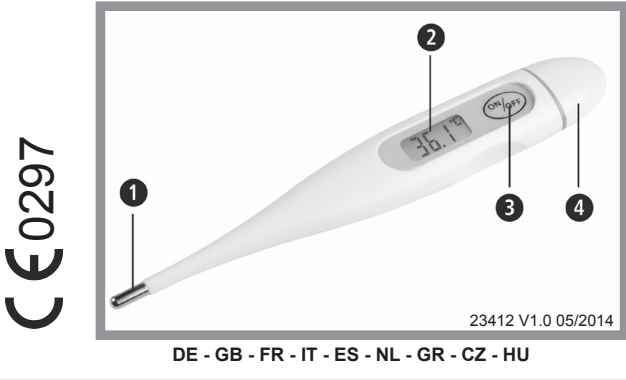## GB Instruction manual Thermometer TM-62E

Thank you for your confidence and congratulations!  
You have purchased a quality product by **ecommed**. To ensure the best results and long-term satisfaction with your **ecommed TM-62E** Thermometer, we recommend that you read the following operating and maintenance instructions carefully.

- 1. Sensor tip for measuring temperature
- 2. LCD screen
- 3. Start button (ON / OFF)
- 4. Battery compartment cover

### SAFETY INFORMATION

- The thermometer contains small parts that could be swallowed by children. Do not let children handle the thermometer without adult supervision.
- Only use the thermometer for measuring body temperature.
- The unit is not designed for commercial use.
- If a fault occurs, please do not try to repair the unit yourself. Discontinue further use of the unit. Contact your dealer and only have repairs carried out by ecommed service personnel.
- Please take note of the operating and storage conditions to ensure that your **ecommed TM-62E** Thermometer will continue to function correctly for a long time.

### GETTING STARTED - Inserting / removing battery:

The **ecommed TM-62E** Thermometer is ready to use right away, as the battery is already installed in the unit. However, if the  symbol should appear on the right of the display after the unit has been in operation for some time, you need to replace the battery. To do this, remove the battery compartment cover **4** and lift the compartment approx. 18 mm out of the thermometer. Use a suitable screwdriver to remove the battery and replace it with a new one of the same type. Make sure the battery is inserted the right way up („+“ - pin is visible). Replace the battery compartment cover.

### OPERATION - Taking a reading:

Press the start button **3** once to switch the unit on. All display symbols will appear for approx. two seconds and you will hear a short tone. After that the last reading taken will be shown shortly in the LCD display **2** (a small „M“ appears on the right side additionally at the same time). If the temperature is below 32 °C, you will see an Lo and a flashing °C on the LCD display **2**. The **ecommed TM-62E** thermometer is now ready to take a temperature reading. Insert the tip of the sensor **1** either orally into the armpit or into the rectum. The thermometer will now automatically begin to take a reading.

During the reading, the current temperature will be shown continuously. As soon as the reading is finished, you will hear 10 beeps and the reading taken will be shown on the display **2**. If the reading is  $> 37,7^{\circ}\text{C}$ , the thermometer will indicate fever alarm by announcing 10 triple beep sounds. The thermometer will switch off automatically eight minutes after taking the reading. To prolong the life of the battery, however, you should switch off the thermometer manually by pressing start button **3** once you have taken the reading.

### Taking temperature readings at various parts of the body

The **ecommed TM-62E** Thermometer is suitable for oral (in the mouth), axillary (in the armpit) and rectal (in the rectum) temperature readings. It is com pletely waterproof. After approx. 1 minute you will hear 10 beep sounds and the reading taken will be shown on the display **2**. You may now remove the thermometer. To ensure safe and accurate measurement, you should take account of the following information when positioning the thermometer:  
ORAL READINGS: Place the sensor tip **1** of the thermometer under the tongue and close your mouth. Breathe through your nose while taking the reading in order to avoid inaccurate results.  
RECTAL READINGS: Many patients find rectal measurement uncomfortable, but it delivers precise results. Carefully insert the sensor tip **1** approx. 3,5 cm into the rectum.  
ARMPIT READINGS: Wash the armpit before taking a reading. Insert the tip of the thermometer **1** into the armpit. Lower the arm and press it lightly against the body while the reading is being taken.

### MEMORY FUNCTION

The **ecommed TM-62E** Thermometer features a memory function that enables you to recall the last reading taken. Press the start button **3** once to switch the unit on. All display symbols will appear for approx. two seconds and you will hear a short tone. After that the last reading taken will be shown shortly in the display **2** (a small „M“ appears on the right side additionally at the same time).

### CLEANING AND DISINFECTION

The **ecommed TM-62E** Thermometer is waterproof, making it easy to clean and disinfect. Clean the thermometer before and after each use using a cloth moistened in a mild solution containing alcohol. You may also soak the thermometer in a suitable disinfectant solution. Your dealer will be happy to recommend suitable cleaning products and give you advice on their use.

### SPECIFICATIONS

**Name and model:** **ecommed TM-62E** Thermometer  
**Measuring range / Scale:** 32,0°C to 42,0 °C / 0,1°C  
**Measuring accuracy:**  $\pm 0,1^{\circ}\text{C}$  (deviation from actual body temperature in the range of 35,5°C to 42°C), non-extrapolating approx. 2.500 readings  
**Operating life:** 10°C to 40°C / -10°C to +60°C  
**Operating / Storage conditions:** 10°C to 40°C / -10°C to +60°C  
**Dimensions / Weight:** 10 x 128 x 19 mm (H x L x W) / approx. 9,8 g  
**Battery:** 1 button cell 1,5 V (type AG3 or LR41)  
**Article No. / EAN Code:** 23412 / 40 15588 23412 5

-  Read instructions carefully before use!
-  LOT number
-  Device complies with protection class B
-  Manufacturer
-  Date of manufacture

The CE-mark refers to the EC-directive 93/42/EEC. The product was designed and manufactured in conformance with the following European norms: EN 12470-3 (Medical thermometers), EN 980 (Symbols), EN 60601-1 (Electrical safety for medical items). **Electromagnetic compatibility:** The device complies with the EN 60601-1-2 standard for electromagnetic compatibility. We reserve the right to make technical and design changes in the course of continuous product improvement.

### DISPOSAL

 This product must not be disposed together with the domestic waste. All users are obliged to hand in all electrical or electronic devices, regardless of whether or not they contain toxic substances, at a municipal or commercial collection point so that they can be disposed of in an environmentally acceptable manner. Please remove the battery before disposing of the equipment. Do not dispose of old batteries with your household waste, but dispose of them at a battery collection station at a recycling site or in a shop. Consult your municipal authority or your dealer for information about disposal.

### WARRANTY

Please contact your dealer or the service centre in case of a claim under the warranty. If you have to send in the unit, please enclose a copy of your receipt and state what the defect is. The following warranty terms apply:

- The warranty period for **ecommed** products is 2 years from date of purchase. In case of a warranty claim, the date of purchase has to be proven by means of the sales receipt or invoice.
- Defects in material or workmanship will be removed free of charge within the warranty period.
- Repairs under warranty do not extend the warranty period either for the unit or for the replacement parts.
- The following is excluded under the warranty:  
**a.** All damage which has arisen due to improper treatment, e.g. nonobservance of the user instructions.  
**b.** All damage which is due to repairs or tampering by the customer or unauthorised third parties.  
**c.** Damage which has arisen during transport from the manufacturer to the consumer or during transport to the service centre.  
**d.** Accessories which are subject to normal wear and tear (batteries etc.).
- Liability for direct or indirect consequential losses caused by the unit is excluded even if the damage to the unit is accepted as a warranty claim.

 **MEDISANA AG, Jagenbergstr. 19, 41468 NEUSS, GERMANY.**  
**ecommed** is a registered trademark of **Medisana AG**.

*For service, accessories and spare parts, please contact:*  
**MEDHEALTH SUPPLIES LTD.**  
Unit 6, Dock offices, Surrey Quays Road, London SE16 2XU  
Tel.: + 44 / 207 - 237 88 99, Fax: + 44 / 207 - 231 13 00  
eMail: info@medhealthcare.co.uk; Internet: www.medisana.co.uk

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